

AXIS 350 **Advanced Line**

Mobile Operating Table with sliding top
for all surgical applications with wide
range of Accessories



Axis 350 – Technical Specifications

Axis 350 - Specifications	Specifications
Length with Head Rest, Long Intermediate Plate and Leg Plates	2185 mm
Length with Head Rest, Short Intermediate Plate and Leg Plates	2035 mm
Width of table top (without side rails)	520 mm
Width of table top (with side rails)	580 mm
Length table pad (without dedicated accessories)	630 mm
Width table pad (without dedicated accessories)	520 mm
Column moving range (high & low)	370 mm
Sliding table top	300 mm
Column Range axis 350 High (lowest & highest position, without pad)	750 - 1120 mm
Column Range axis 350 Low (lowest & highest position, without pad)	612 - 982 mm
Pad thickness	69 mm
X-ray window A-Position	1285 x 348 mm
X-ray window B-Position	1135 x 348 mm
Side rail dimension	25 x 10 mm
Table weight (without dedicated accessories)	255 kg
Static workload	450 kg
Dynamic workload	220 kg
Wheels diameter	125 mm
Seat Plate up / down	80° / 40°
Trendelenburg / Reversed Trendelenburg	30°
Lateral tilting left / right	20°
Head Rest up / down manual	45°
Leg Plate up / down manual	-90° / 20°
IPX – Protection against entry of liquids	IPX4

Batteries

➤ Main Power Supply

- ✓ 2 Rechargeable Lead Acid Batteries
- ✓ 22 Ah, 12V (24V as System)
- ✓ Lifetime: 500 Cycles approx. 2 Years
- ✓ 50 Operations, capacity of appr. One week usage

➤ Emergency Power Supply

- ✓ 2 Rechargeable Lead Acid Batteries
- ✓ 3.2 Ah, 12V (24V as System)
- ✓ Lifetime: 500 Cycles approx. 2 Years










➤ MOP Power Supply

- ✓ Rechargeable Lithium-Ionen Battery
- ✓ 2.6 Ah, 3.6V
- ✓ Lifetime: 500 Cycles approx. 2 Years



❖ For long battery life – charge table every day over night!


Axis 350 - Dimensions

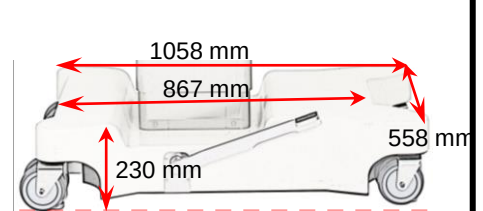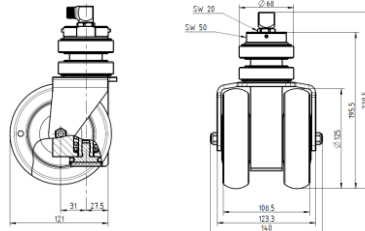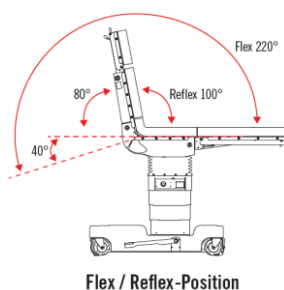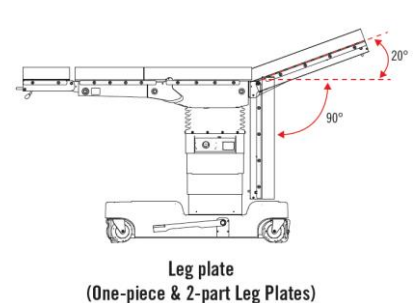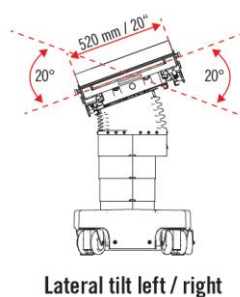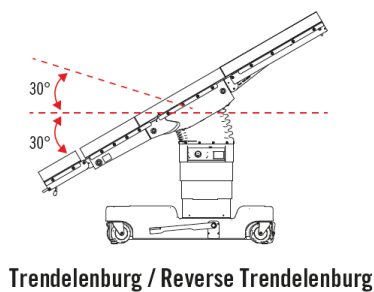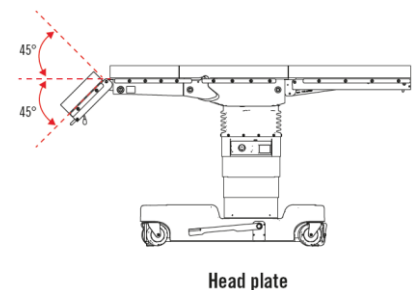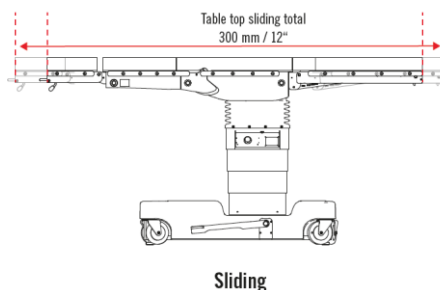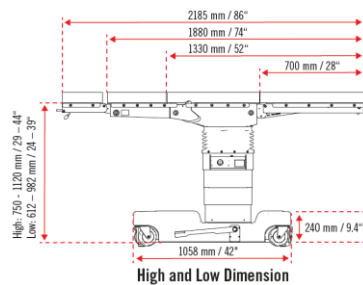

Foot dimensions


Components - Base

- 1 - Head plate
- 2 - Intermediate plate
- 3 - Leg plate
- 4 - Emergency operational panel
- 5 - Socket for wired hand control
- 6 - Lever for foot pump
- 7 - Terminal for potential equalization
- 8 - Socket for control device
- 9 - Socket for main cable
- 10 - Battery charging indicator

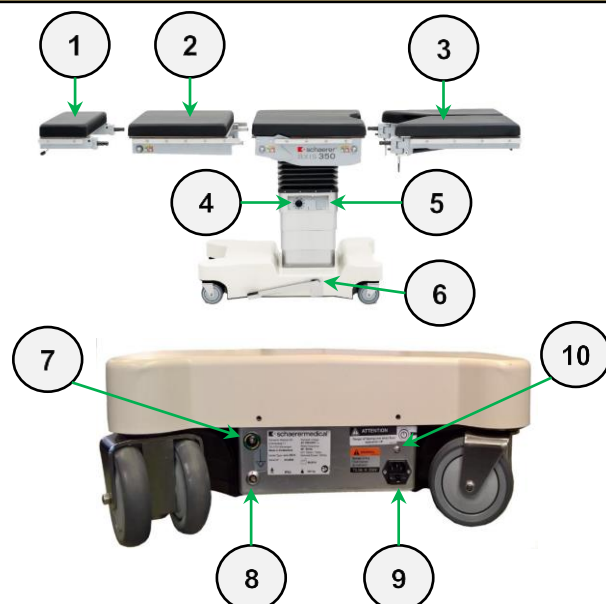

Axis 350 – Technical Specifications

Axis 350 – Segments

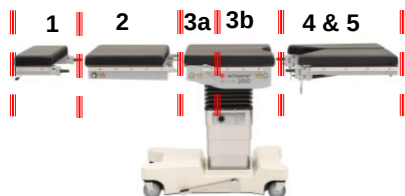

Axis 350 – Sliding Top

Manual sliding table top **300 mm**
with **electric** lock/unlock button



Axis 350 - Range


Low Version

- Without pads: 612 to 982 mm
- With pads: 681 to 1051 mm
- Height pads: 69 mm



High Version

- Without pads: 750 to 1120 mm
- With pads: 819 to 1189 mm
- Height pads: 69 mm

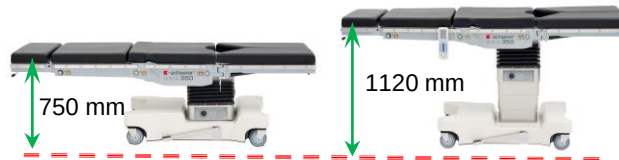

Axis 350 – Length

Head Rest, **Long** Intermediate Plate and Leg Plates



Head Rest, **Short** Intermediate Plate and Leg Plates



Axis 350 – Patient Weight Capacity

Static

450 kg
Static



Dynamic **A-Position** full slide

185 kg
Dynamic



Dynamic **B-Position** full slide

220 kg
Dynamic



Axis 350 – X-Ray window A-Position



Axis 350 – X-Ray window B-Position



Axis 350 – Hand Control

- Wired Hand Control (standard) with two connecting points on the operating table
- or
- Wireless Hand Control (optional) with RFID and Bluetooth technology

➤ Features:

- ✓ Colored backlit touch display 71 mm (2.8 inch)
- ✓ One Hand Control for multiple tables
- ✓ Individually programmable memory function
- ✓ Error messages and Log-file for service purposes
- ✓ Liquid protection: IPX7
- ✓ No need of extra charging station
- ✓ Motor drive control. No need for extra drive control unit



Item: **92523**

Wireless Hand control
(optional)

Axis 350 – Data Log on SD Card

Data Log supports failure analysis
Quick response on the failures



Axis 350 – Emergency Control

Manual emergency control system with foot pedal



Axis 350 – Foot Control

Possible Movements:

- ✓ Column up and down
- ✓ Trendelenburg / Reverse Trendelenburg
- ✓ Tilt right and left



Item: **92532**

Wired foot control
Axis 350



Item: **92533**

Wireless foot control
Axis 350

Axis 350 – Components

Axis 350 – High Version

Item: **92600H**
Axis 350
High Version




Axis 350 – Low Version

Item: **92600L**
Axis 350
Low Version




5th Wheel (optional)

Item: **92261**
Direction Control
Axis 350

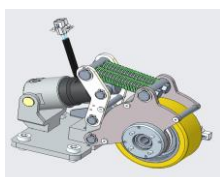


Item: **92260**
Drive Unit
Axis 350

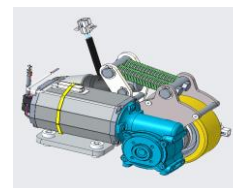

Head Rest

Item: **92640**
Head rest

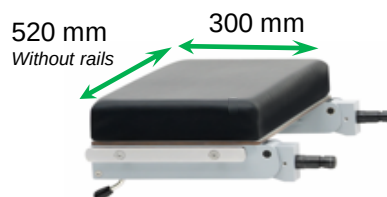


Intermediate Plates


Item: **92660**
Short seat plate



Item: **92631**
Long Intermediate plate

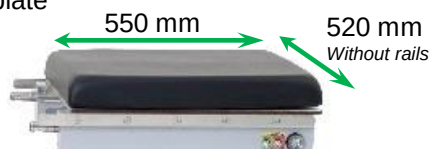


Item: **92630**
Short Intermediate plate

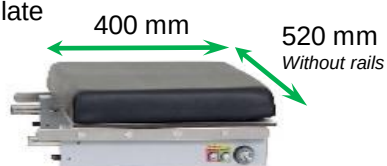


Leg Plates


Item: **92632**
Leg plates 2-part

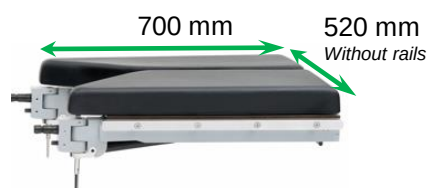

Item: **92633**
Leg plate one-piece

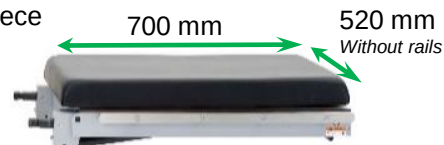

Item: **92620**
Leg plate light

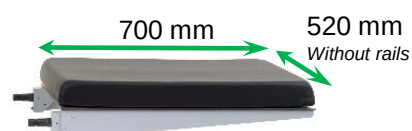

Axis 350 – Components

Connection Combinations




A - Position



B - Position

Combination	Connections						Length			
Combination 1	Head Plate 	+	Int. Plate Long 	+	Axis 350 	+	Leg Plates 	= 2180 mm		
Combination 2	Head Plate 	+	Int. Plate Short 	+	Axis 350 	+	Leg Plates 	= 2030 mm		
Combination 3	Leg Plates 		+	Axis 350 	+	Int. Plate Long 	+	Head Plate 	= 2180 mm	
Combination 4	Leg Plates 		+	Axis 350 	+	Int. Plate Short 	+	Head Plate 	= 2030 mm	
Combination 5	Head Plate 	+	Short Seat 	+	Axis 350 	+	Int. Plate Long 	+	Leg Plates 	= 2359 mm
Combination 6	Head Plate 	+	Short Seat 	+	Axis 350 	+	Int. Plate Short 	+	Leg Plates 	= 2209 mm

Basic Configuration

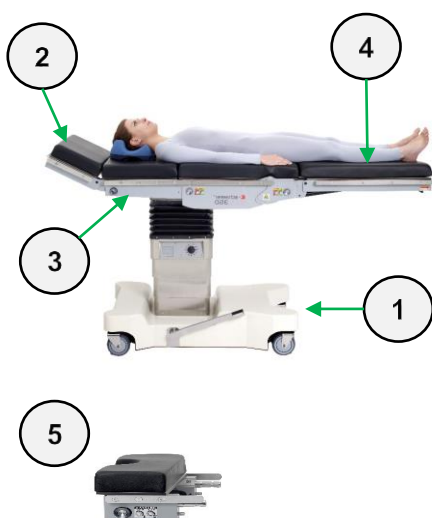

No	Item	Qty	Description
1	92600L 92600H	1	Axis 350 - Low Version (or) Axis 350 - High Version
2	92640	1	Head rest Axis 350
3	92630 92631	1	Short intermediate plate Axis 350 (or) Long intermediate plate Axis 350
4	92632 92633 92620	1	2-part Leg plates Axis 350 (or) Leg plate one-piece (or) Leg plate light
	91260	1	Packaging OT-Table WITHOUT ramp
	91261	1	Reusable exit ramp
Optional Accessories			
5	92660	1	Short seat plate
	92261	1	5th Wheel direction Control

Accessories for General Surgery

6



Item: **90321**

(or)



Item: **90323**

(or)



Item: **90635**

(or)



Item: **90636**

(+)



Item: **90506**

(or)



Item: **90310**

7



Item: **90520**

(or)



Item: **4646.05**
Item: **5.890-VE**

8



Item: **800-0057**

9



Item: **92021**

10



Item: **90500**

No	Item	Qty	Description
6	90321	1	Anaesthesia screen, rigid (+)
	90506	1	Radial setting clamp with fixing screw (or)
	90310		Rapid-action radial setting clamp
6	90323	1	Anaesthesia screen with central articulation (+)
	90506	1	Radial setting clamp with fixing screw (or)
	90310		Rapid-action radial setting clamp
6	90635	1	Anaesthesia screen, with extendable rods (+)
	90506	1	Radial setting clamp with fixing screw (or)
	90310		Rapid-action radial setting clamp
6	90636	1	Anaesthesia screen central art., ext. rods (+)
	90506	1	Radial setting clamp with fixing screw (or)
	90310		Rapid-action radial setting clamp
7	90520	2	Arm rest (or)
	4646.05		Arm support quick central fixation (+)
	5.890-VE		Fastening straps (1 pair)
8	800-0057	1	Restraint Strap
9	92021	1	Arm strap with fixing screw
10	90500	1	Infusion pole with clamp, height-adjustable
11	90332	1	Instrument tray with clamp, height-adjustable
12	90360	1	Accessories trolley (+)
	90300	1	Packaging for Accessories trolley
Dedicated Accessories			
13	92029	1	Universal Adapter for DORO® basic unit
14	90939	1	Table extension
15	90913	1	Hand surgery table with fixing screw
16	92006	1	X-ray cassette frame for Axis 350


Ophthalmology / ENT – Configuration A

1



Item: 92029

2

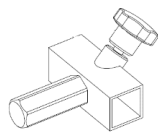

Item: 90456

3




Item: 90041

4



Item: 90003

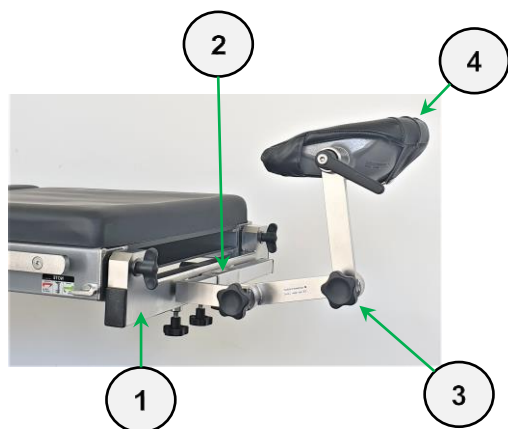

No	Item	Qty	Description
1	92029	1	Adapter for DORO® basic unit, for Axis 350
2	90456	1	Hexagonal adapter with square adapter
3	90041	1	Head rest holder, with lateral adjustment
4	90003	1	Head rest (horseshoe) with pad

Ophthalmology / ENT – Configuration B

1



Item: 92029

2



Item: 90045

3



Item: 1002.72A0


4



Item: 1004.75A0

No	Item	Qty	Description
1	92029	1	Adapter for DORO® basic unit, for Axis 350
2	90045	1	Horseshoe-headrest holder 3 joints
3	1002.72A0	1	Horseshoe-shaped head rest, two pieces
4	1004.75A0	1	Wrist hand support for ophthalmology

Abdominal – Configuration A

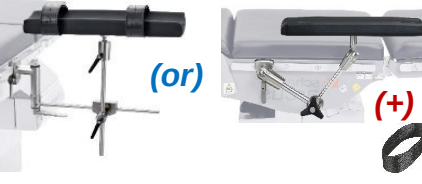
1  (or) (+) Item: 90520 Item: 4646.05 Item: 5.890-VE	2  (+) (or) Item: 90325 Item: 90506 Item: 90310	3  (+) (or) Item: 90330 Item: 90506 Item: 90310
---	---	--



No	Item	Qty	Description
1	90520	1	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
2	90325 90506	1 1	Arm strap with fixing screw Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp
3	90330 90506	2 2	Leg holder "Goepel" (1 piece) (+) Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp


Abdominal – Configuration B

1  (or) (+) Item: 90520 Item: 4646.05 Item: 5.890-VE	2  (+) (or) Item: 90325 Item: 90506 Item: 90310	3  (+) Adjustable Leg holder
---	---	--

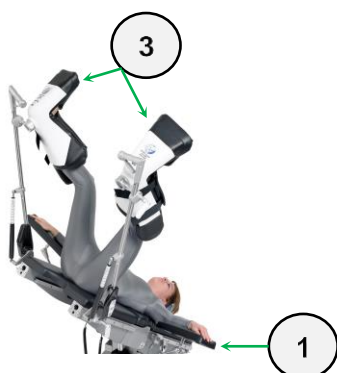

No	Item	Qty	Description
1	90520	1	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
2	90325 90506	1 1	Arm strap with fixing screw Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp
3	800-0342-PR 800-0338-EU	1 2	Adjustable Leg holder for patients up to 182 kg (or) 800-0342-PL for patients up to 272 kg (or) 800-0364 for patients up to 454 kg (+) Springloc Clamp (order 2 pieces)

Bariatric, Laparoscopic – Configuration A

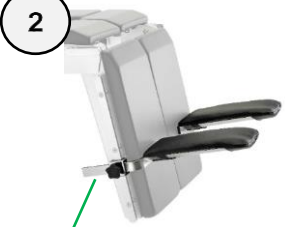1  Item: 90520 (or)  Item: 4646.05 + 5.890-VE	2  Item: 90079 (+)  Item: 81501	3  Item: 800-0055 (+)  Item: 800-0338-EU	4  Item: 90939
--	--	--	---



No	Item	Qty	Description
1	90520	1	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
2	90079 81501	1 2	Knee cradles (1 pair), with pads (+) Clamp, Knee cradles (order 2 pieces)
Optional			
3	800-0055 800-0338-EU	1 2	Nissen straps, 272 kg patient w. cap. (1 pair) (+) Springloc Clamp
4	90939	4	Table extension (1 piece),

Bariatric, Laparoscopic – Configuration B

1  Item: 90520  Item: 4646.05 + 5.890-VE	2  (+)  Adjustable Leg holder	3  Item: 800-0055 (+)  Item: 800-0338-EU	4  Item: 90939
---	--	--	---



No	Item	Qty	Description
1	90520	1	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
2	800-0342-PR 800-0338-EU	1 2	Adjustable Leg holder for patients up to 182 kg (or) 800-0342-PL for patients up to 272 kg (or) 800-0364 for patients up to 454 kg (+) Springloc Clamp (order 2 pieces)
3	800-0055 800-0338-EU	1 2	Nissen straps, 272 kg patient w. cap. (1 pair) (+) Springloc Clamp
4	90939	4	Table extension (1 piece),

Nephrectomy, Thorax

1



Item: 90520

(or)



Item: 4646.05
Item: 5.890-VE

(+)



2



Item: 90325

(+)



Item: 90506

(or)



Item: 90310

3



Item: 90351

(+)



Item: 90350

(+)



Item: 90539

(or)



Item: 90314

4



Item: 90352

(+)



Item: 90350

(+)



Item: 90539

(or)



Item: 90314

5



Item: 90609

(+)



Item: 90350

(+)




Item: 90539

(or)



Item: 90314

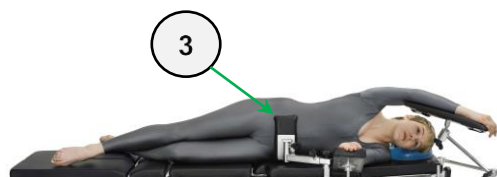

No	Item	Qty	Description
1	90520	2	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
2	90325 90506	1 1	Arm rest high, without clamp (+) Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp
3	90351 90350 90539	1 1 1	Side and back rest, with pad (+) Universal holder for side and back rest (+) Radial clamp with fixing screw, rotating (or) 90314 Rapid-action radial setting clamp
4	90352 90350 90539	1 1 1	Side and back supporting roll Universal holder for side and back rest (+) Radial clamp with fixing screw, rotating (or) 90314 Rapid-action radial setting clamp
5	90609 90350 90539	1 1 1	Sacrum/Pubis/Sternum support, with pad (+) Universal holder for side and back rest (+) Radial clamp with fixing screw, rotating (or) 90314 Rapid-action radial setting clamp

Rectal, Andrews, Genupectoral

1



Item: 90520

(or)



Item: 4646.05
Item: 5.890-VE

2



Item: 90277

(+)



Item: 90506

(or)


Item: 90310

3



Item: 90079

(+)

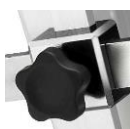

Item: 81501

4




Item: AP022

(or)



Item: QPHPD-ADULT 1

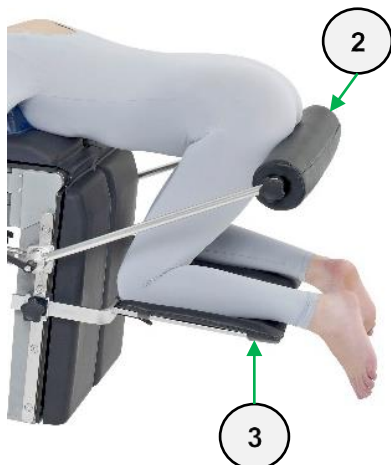



No	Item	Qty	Description
1	90520	2	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
2	90277 90506	1 2	Abdominal rest, with pad, adjustable (+) Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp
3	90079 81501	1 2	Knee cradles (1 pair), with pads (+) Clamp, Knee cradles (order 2 pieces)
4	AP022	1	Gel Headrest for prone position (or) QPHPD-ADULT 1 – Head rest (single use)

Gynaecology / Urology – Configuration A




No	Item	Qty	Description
1	90330 90506	2 2	Leg holder "Goepel" (1 piece) (+) Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp
2	9911202 9911012	1 1	Tray TUR (+) Frame for Tray TUR
3	90520	2	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)


Gynaecology / Urology – Configuration B

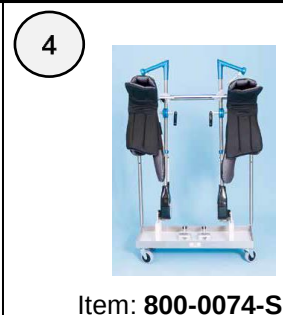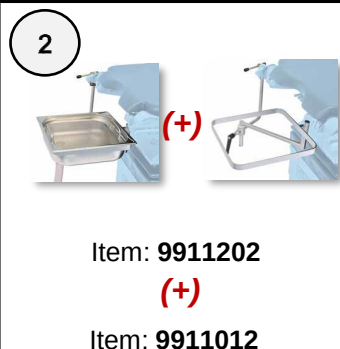

No	Item	Qty	Description
1	800-0342-PR 800-0338-EU	1 2	Adjustable Leg holder for patients up to 182 kg (or) 800-0342-PL for patients up to 272 kg (or) 800-0364 for patients up to 454 kg (+) Springloc Clamp (order 2 pieces)
2	9911202 9911012	1 1	Tray TUR (+) Frame for Tray TUR
3	90520	2	Arm rest with clamp (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
Optional			
4	800-0074-S	1	Great White® Stirrup Dolly



Spine – Configuration

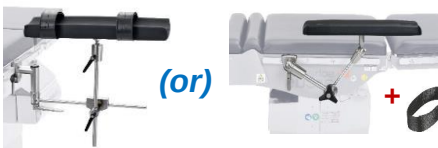1  Item: NO 165	2  Item: NA 401	3  Item: 90520 Item: 4646.05 Item: 5.890-VE
4  Item: 92029 Item: 90456 Item: 90041 Item: 90003 Item: QPHPD-ADULT 1 Item: AP022		
5  Item: 800-0235	6  Item: 800-0074-SF	

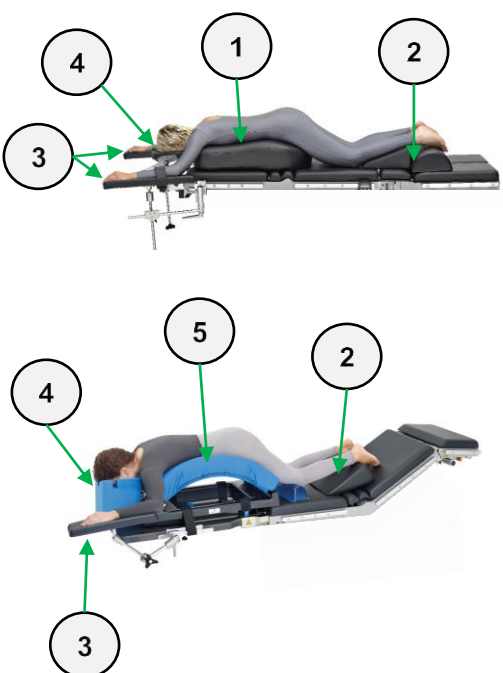


No	Item	Qty	Description
1	NO 165	1	Torso Pad "Onsernone"
2	NA 401	1	Leg support wave cushion 65x50x13cm
3	90520	2	Arm rest (or) 4646.05 Arm support quick central fixation (+) 5.890-VE Fastening straps (1 pair)
4	92029 90456 90041 90003	1 1 1 1	Adapter for DORO® basic unit, for Axis 350 (+) Hexagonal adapter with square adapter (+) Head rest holder, with lateral adjustment (+) Head rest (horseshoe) with pad (or) QPHPD-ADULT 1 - QView™ Head positioning device (single use) (or) AP022 Gel Headrest for prone position
Optional Accessories			
5	800-0235	1	SchureSpine Frame, Radiolucent
6	800-0074-SF	1	Schure Spine Frame Dolly

Neuro – Configuration A - Aluminium Skull Clamp System

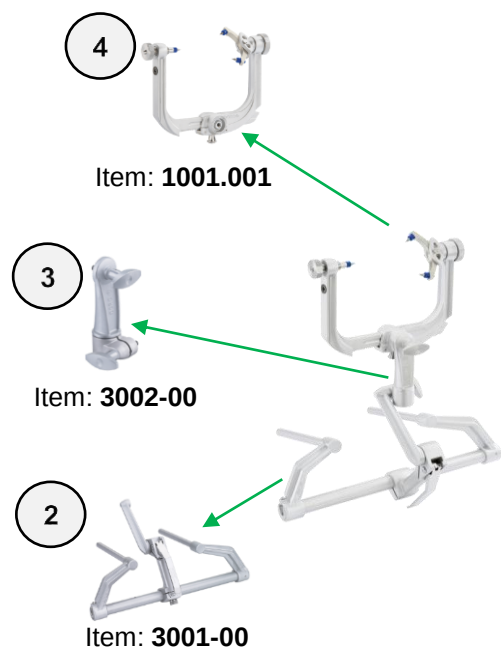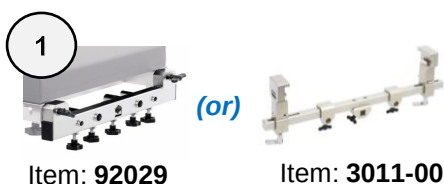

No	Item	Qty	Description
1	92029	1	Adapter for DORO® basic unit, Axis 350 (or) 301110 Rail adapter (for tables not Schaerer)
2	3001-00	1	DORO® Base unit aluminum
3	3002-00	1	DORO® Swivel adaptor aluminum (or) DORO® Ball adapter, allows 360° rotation
4	1001.001	1	DORO® Skull clamp QR3 aluminum
5	3009-00	1	DORO® Horseshoe headrest for adults
6	3008-00	1	DORO® Horseshoe headrest for children
		1	Pins reusable (or) Pins disposable
Additional accessories for sitting position			
7	3007-00 90506	1	DORO® Rail adapter for sitting position (+) Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp

Sitting Position

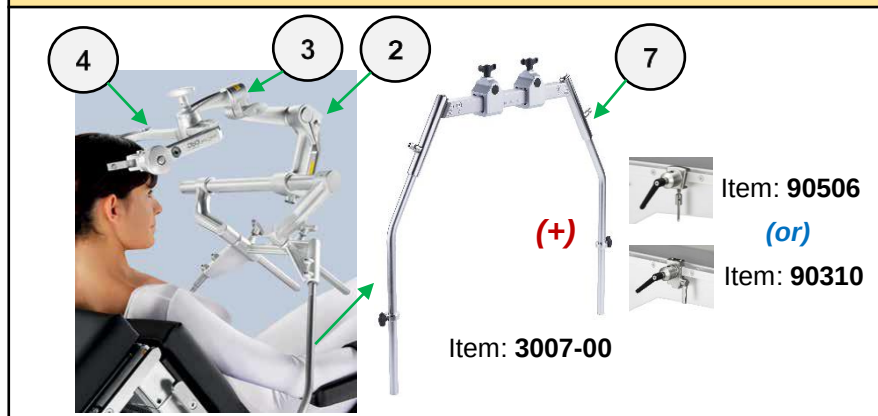

Pins reusable



Item: **3005-00**
Reusable skull pins
for **adults** (3 pieces)



Item: **3004-00**
Reusable skull pins for
children (3 pieces)

Pins disposable



Item: **1106.003**
DORO® disposable skull
pins, for **adults**,
(36 pieces)



Item: **3006-10**
DORO® disposable
skull pins, for **children**,
(36 pieces)

Pins disposable



Item: **3006-20**
DORO® disposable
skull pins, for **adults**,
(36 pieces)





Item: **3006-30**
DORO® disposable
skull pins, for **children**,
(36 pieces)

Neuro – Configuration B - Radiolucent Skull Clamp System


Item: **1101.020** (Complete Set)

DORO® Radiolucent Head Rest System compatible with all tables
(No adapter needed)

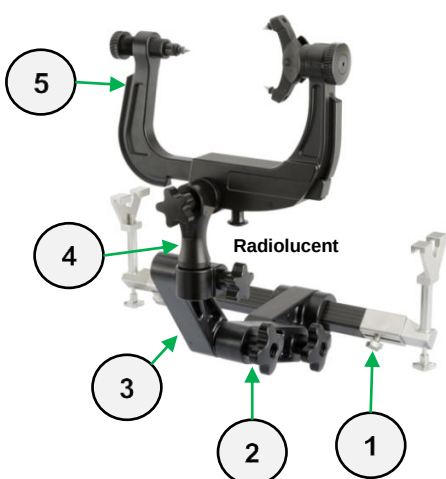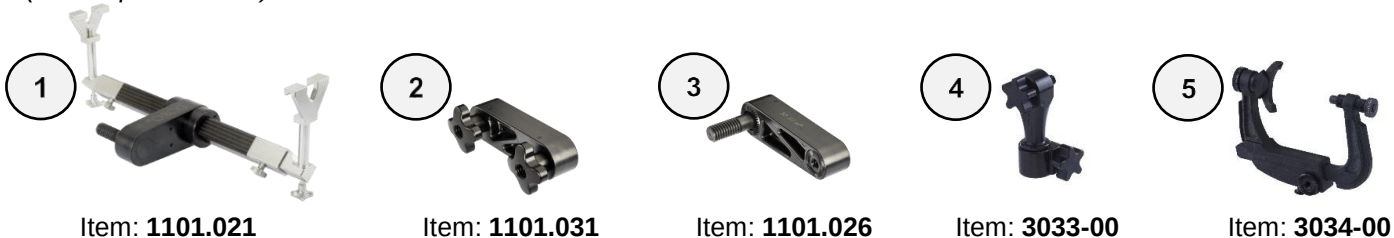

No	Item	Qty	Description
1	1101.021	1	DORO® LUCENT Basic support
2	1101.031	1	DORO® LUCENT Clamp intermediate piece
3	1101.026	1	DORO® LUCENT Intermediate piece
4	3033-00	1	DORO® Swivel adapter, radiolucent
5	3034-00	1	DORO® Skull clamp, radi. with Quick-Rail™
6		1	Pins reusable (or) Pins disposable
Additional accessories for sitting position			
7	3007-00 90506	1	DORO® Rail adapter for sitting position (+) Radial setting clamp with fixing screw (or) 90310 Rapid-action radial setting clamp
8	3001-00	1	DORO® Base unit aluminum

Sitting Position

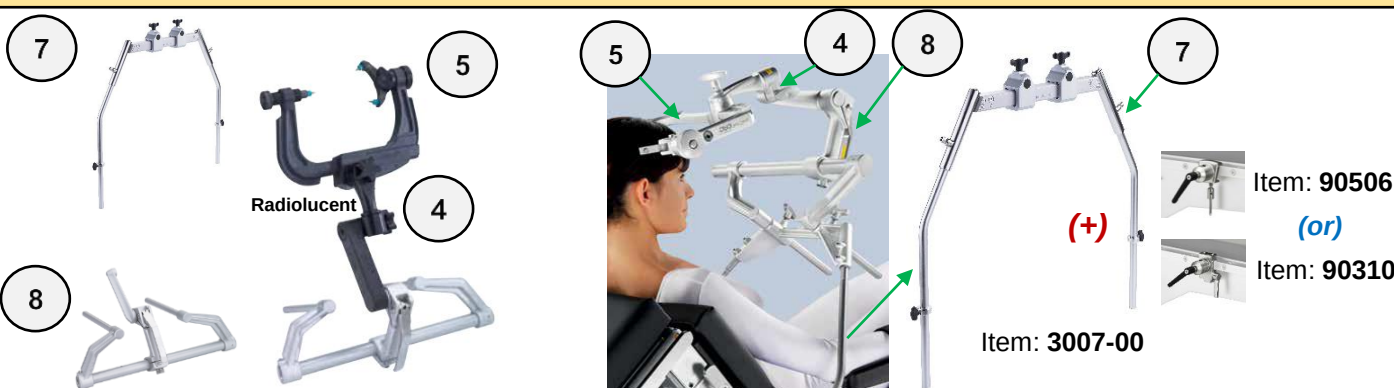

Pins reusable



Item: **3005-00**
Reusable skull pins
for **adults** (3 pieces)



Item: **3004-00**
Reusable skull pins for
children (3 pieces)

Pins disposable



Item: **1106.003**
DORO® disposable skull
pins, for **adults**,
(36 pieces)



Item: **3006-10**
DORO® disposable
skull pins, for **children**,
(36 pieces)

Pins disposable



Item: **3006-20**
DORO® disposable
skull pins, for **adults**,
(36 pieces)



Item: **3006-30**
DORO® disposable
skull pins, for **children**,
(36 pieces)

ETS – Traction System

Applications

Femur






Tibia



Lateral Position




ETS - Femur

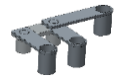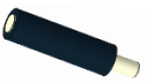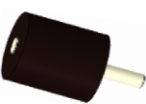No	Images	Item	Description	No	Images	Item	Description
1		92690	Connection bar ETS	6		92693	Bar ETS (Changeable from left to right - order 2 pc.)
2		92691	Seat plate ETS	7		92694	Traction Unit ETS (Changeable from left to right - order 2 pc.)
3		90958	Offset holder for Perineal rod	8		92695	ETS Right Leg plate
4		90971	Perineal rod	9		92696	ETS Left Leg plate
5		90990	Perineal rod for Arthroscopy	10	 <i>or</i> 	90763 2031113 2031470	Extension shoe white (Changeable from left to right - order 2 pc.) <i>or</i> Footplate ZSB, right + Footplate ZSB, left

ETS - Tibia Additional Items

No	Images	Item	Description
11		92612	Traction Unit for Tibia ETS
12		92699	Foot plate
13		92698	Foot Plate Child
14		90193	Telescopic support

Optional Items

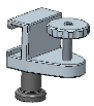No	Images	Item	Description
16		90520	Arm rest with rapid-action fixing clamp
17		90325 90506	Arm rest high (+) Radial setting clamp
18		90330 90506	Goepel (+) Radial setting clamp
19		81696	Side rail extension
20		90994	Extension bow holder


ETS - Lateral Position Additional Items

No	Images	Item	Description
15		92611	ETS Counter Traction Unit for lateral position

ETS – Trolley


No	Images	Item	Description
21		92613	ETS Trolley

ETS - 1 Leg Traction

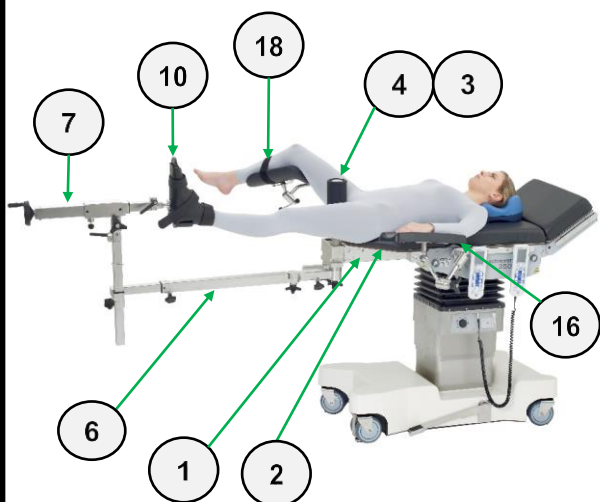


No	Item	Qty	Description
1	92690	1	Connection bar ETS
2	92691	1	Seat plate ETS
3	90958	1	Offset holder for perineal rod (optional)
4	90971	1	Perineal rod
6	92693	1	Bar ETS
7	92694	1	Traction Unit ETS
8	92695	1	ETS Right Leg Plate
9	92696	1	ETS Left Leg Plate
10	90763 2031113 2031470	1	Extension shoe white (or) Footplate ZSB, right (+) Footplate ZSB, left
16	90520	2	Arm rest with rapid-action clamp
18	90330 90506	1	Goepel (+) Radial setting clamp
20	90994	1	Extension bow holder (optional)

ETS - 2 Leg Traction

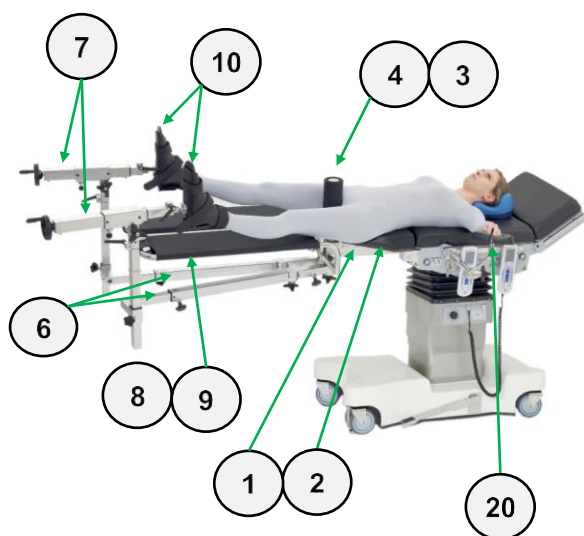


No	Item	Qty	Description
1	92690	1	Connection bar ETS
2	92691	1	Seat plate ETS
3	90958	1	Offset holder for perineal rod (optional)
4	90971	1	Perineal rod
6	92693	2	Bar ETS
7	92694	2	Traction Unit ETS
8	92695	1	ETS Right Leg Plate
9	92696	1	ETS Left Leg Plate
10	90763 2031113 2031470	2 1 1	Extension shoe white (or) Footplate ZSB, right (+) Footplate ZSB, left
16	90520	1	Arm rest with rapid-action clamp
20	90994	1	Extension bow holder (optional)

ETS – Tibia (in addition to Femur Traction)

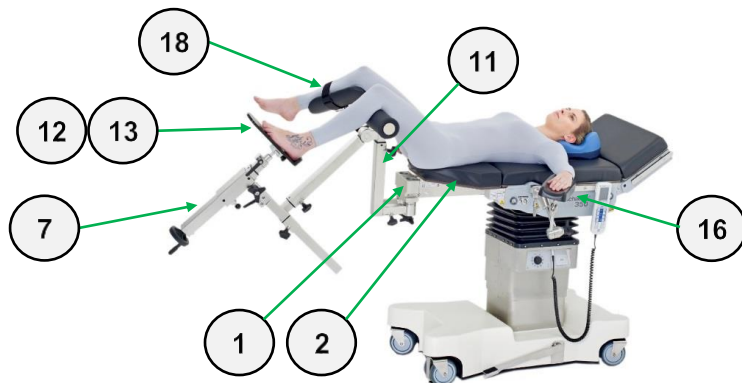

No	Item	Qty	Description
11	92612	1	Traction Unit for Tibia ETS
12	92699	1	Foot plate
13	92698	1	Foot plate Child
14	90193	1	Telescopic support
18	90330 90506	1	Goepel (+) Radial setting clamp
Optional Accessories			
16	90520	2	Arm rest with rapid-action fixing clamp

ETS – Lateral Positioning (in addition to Femur Traction)

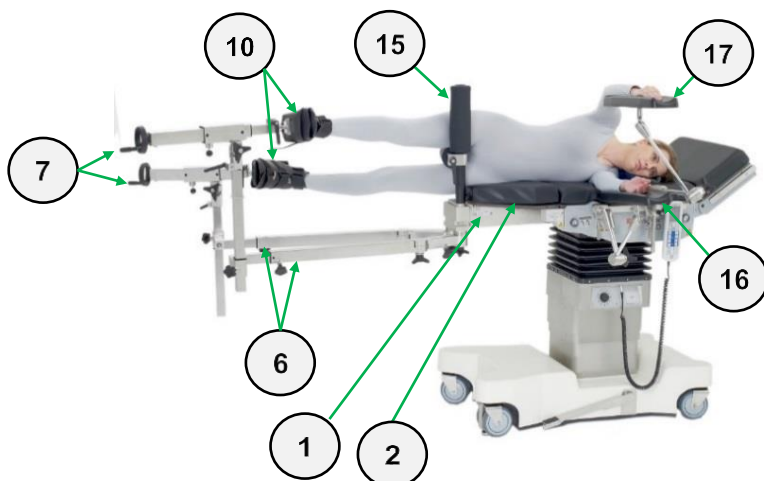

No	Item	Qty	Description
15	92611	1	ETS Counter Traction Unit for lateral position
Optional Accessories			
16	90520	1	Arm rest with rapid-action fixing clamp
17	90325 90506	1 1	Arm rest high (+) Radial setting clamp

Hand Surgery



1



1

(or)



2

(+)



Item: **90506**
(order 2 pc.)

No	Item	Qty	Description
1	90913	1	Hand surgery table, with fixing screw (or)
2	800-0040 90506	1 2	Hand surgery table, radiolucent Radial setting clamp

Humerus / Traction Device



1



No	Item	Qty	Description
1	9907038	1	Humerus and Forearm Positioning Device

Knee Arthroscopy



1



2

(or)



Item: **90539**

(+)



Item: **90314**

(or)



Item: **90506**

(+)



Item: **90310**




1



2

No	Item	Qty	Description
1	90357	1	Leg holder arthroscopy (or)
	90539	1	Knee rest (+)
	90539	1	Radial clamp with fixing screw, rotating (or)
	90314	1	Rapid-action radial clamp, rotating
2	90330	1	Leg holder "Goepel" (1 piece) (+)
	90506	1	Radial setting clamp with fixing screw (or)
	90310	1	Rapid-action radial setting clamp


SSP - Schaerer Shoulder Plate - Configuration A

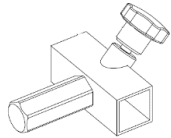
1  Item: 92625	2  Item: 90456	3  Item: 90041	4  Item: 90566 (or) Item: 90003
---	---	--	--



No	Item	Qty	Description
1	92625	1	Shoulder plate (SSP) Axis 350
2	90456	1	Adapter
3	90041	1	Head rest holder, with lateral adjustment
4	90566 90003	1	Helmet, radiolucent (or) Head rest (horseshoe) with pad

SSP - Schaerer Shoulder Plate - Configuration B

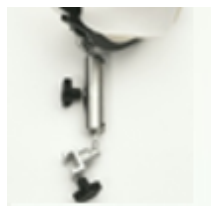1  Item: 92625	2  Item: 92666	3  Item: 90899	4  Item: F-90023
---	---	--	---

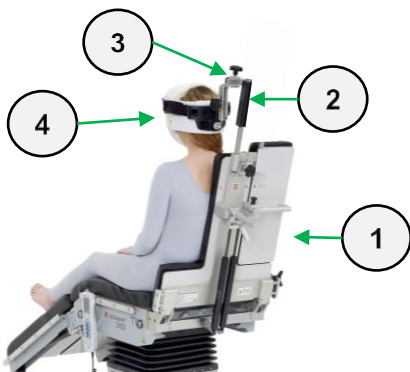


No	Item	Qty	Description
1	92625	1	Shoulder plate (SSP) Axis 350
2	92666	1	Head rest holder
3	90899	1	Head rest
4	F-90023	1	Head pad (10 Pieces)

Axis 350 – One Table – Multiple Features



Data Log
on SD Card

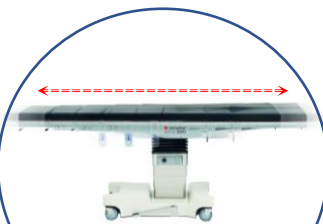

Sliding 300mm



X-Ray Window
A-position: 1285mm



450 kg

Static Weight



Bluetooth wireless
Hand control unit



Emergency Panel



Pairing panel for
wireless
Hand control




Footpump



Quick click buttons






Double wheels 125mm
Cylinder Stamps

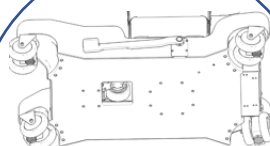

5th wheel
drive unit (optional)



Rubber Bellows

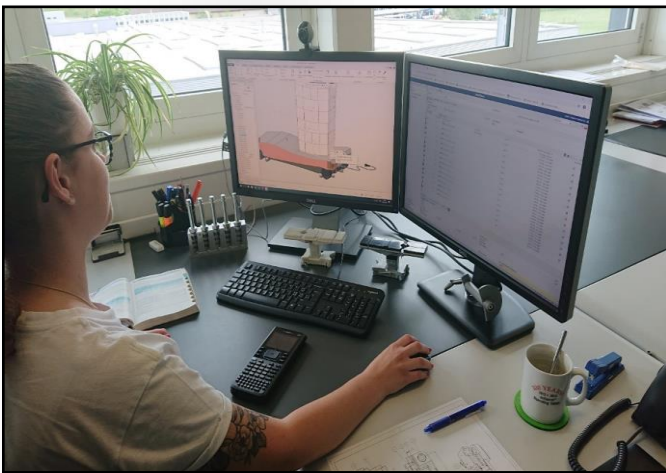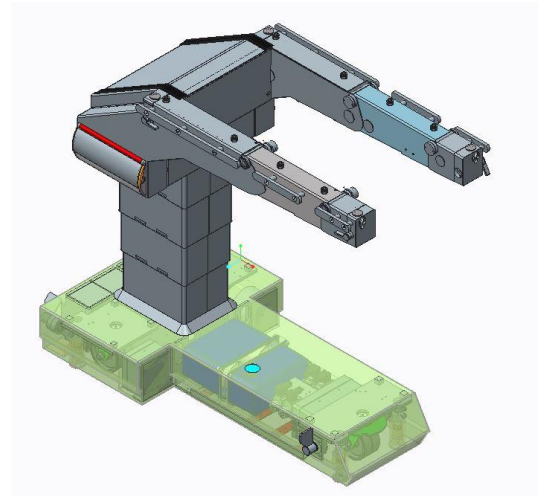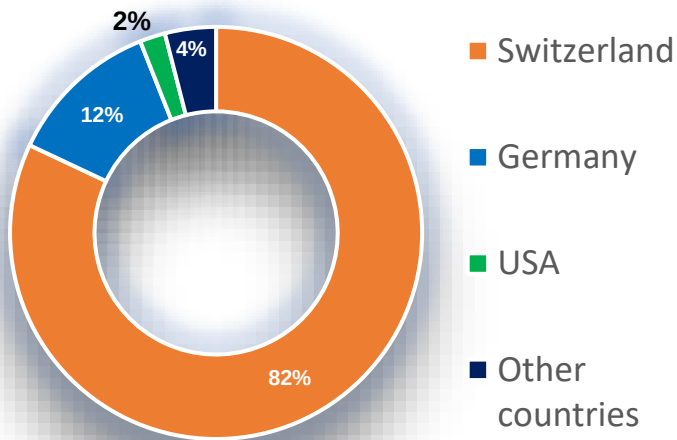

Head office Switzerland



Branch office USA



- ✓ **Design innovation**
Copied often – never reached
- ✓ **Material innovation**
Carbon accessories
- ✓ **Engineering innovation**
Best-in-class features and versatility
- ✓ **Ergonomic innovation**
Tables and accessories designed for ease of use
- ✓ **Safety and hygiene innovation**
Patient safety and hygiene features

Partner network in over 75 countries



